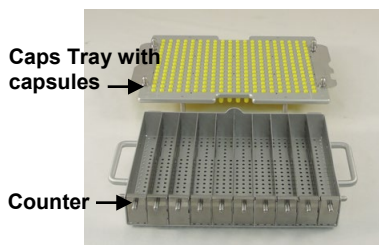




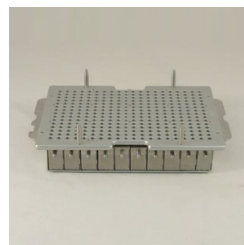
These innovative capsule counters divide a batch of 300 filled capsules directly from the capsule filler's Caps Tray into 10 lots of 30 capsules; 5 lots of 60 capsules; or 3 lots of 100 capsules.



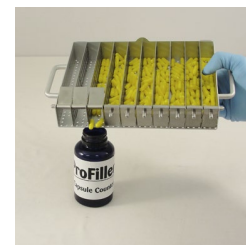
1. After filling & locking capsules, remove Caps Tray from Filler.



2. Place Counter open side down on Caps Tray. Hold them together and turn over.


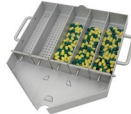
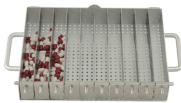


3. Capsules are now divided into smaller lots. Remove Caps Tray



4. Open individual gates (push to left) and tilt Counter to transfer capsules into a bottle.

- Compatible with ProFiller (3600, 3700 & 3800) & ProFill (3006, 3007 & 3008) 300 hole fillers.
- Does not fit other brands of capsule filling machines.
- For bottles with less than 30mm (1 1/8") opening, a funnel may be required for easier transfer.
- Made from anodized aluminum.
- **Tip:** Counter is perforated to allow loose powder to fall through. Shake counter prior to filling bottles to clean capsules of loose powder.

			
Item #	PX CTR 3X100	PX CTR 5X60	PX CTR 10X30
Weight gm (oz)	550 (20)	750 (27)	800 (28)
Size LxWxH cms (inches)	42x18x4 (17x7x1.5)	31x29x4 (12x11x1.5)	31x19x4 (12x7.5x1.5)

